

CELL-LAV

Sterile Tubing Welder



Sterile Tubing Welder


CL30-1S



<< Specially designed for
cellular & genetic drug enterprise partners

Product Overview

CELL-LAV Serile Tubing Welder is widely used for blood filling and preparation in blood stations, hospital blood banks, as well as tubing connection and liquid filling in bioprocessing companies and immunotherapy centers. It provides different welding temperature settings for various bag sizes for different purposes, to meet the tubing welding requirements of different materials such as PVC, latex, silicone and PE tubing.

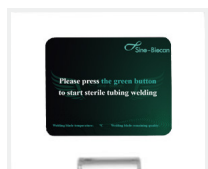

①

Real-time display welding blade temperature
Real-time display remaining welding blade amount



②

Welding blade box
100 pieces/box



③

Applicable to tubing with
Outer diameter 3.9mm–6.0mm
Inner diameter 2.9mm–4.1mm



Applications

- Aseptic connection of 2 blood bag tubings.
- Medical standard tubings with outer diameter 3.9mm–6.0mm and inner diameter 2.9mm–4.1mm.
- Medical PVC tubing.
- Widely applied in: component mixing, leukocyte removal, pathogen inactivation, component dispensing, quality control sampling, cell washing and freezing.

Technical Advantages

- Easy to operate, compact and portable.
- Independent welding blade as consumables, 100 pieces per box; real-time display welding blade temperature and remaining amount, welding blade direct heating simplifies blade production process.
- Connect tubings aseptically in a common non-sterile environment.
- Connect tubings of the same or different tubing diameters, applicable in dry-dry, dry-wet, wet-wet tubing connections.
- The composition and shelf life of the contents in blood bags connected by the aseptic tubing welder are not changed.

Sino-Biocan (Shanghai) Biotech Ltd

Tel:021-60890755 / 4000056700

E-mail:info@sinobiocan.com

Website:www.sinobiocan.net

Address: 2nd Floor, Building 15, No. 6055, Jinhai Road, Fengxian District, Shanghai

Sino-Biocan Biotech (Shanghai) Ltd

Tel:021-60890755 / 4000056700

E-mail:info@sinobiocan.com

Website:www.sinobiocan.net

Address: 1st Floor, Building 15, No. 6055, Jinhai Road, Fengxian District, Shanghai

Sino-Biocan (Chengdu) Biotech Ltd

Tel:028-85840406 / 85893618

E-mail:info@sinobiocan.com

Website:www.sinobiocan.net

Address:4th Floor, Building 4, Tianfu International Biology City, No.618, Fenghuang Road, Shuangliu District, Chengdu